


# Non-Sterile Sales Syringe Adapter Instructions for Use

## Important

- Do not attempt to use the Sales Syringe Adapter (SSA) before completely reading and understanding the information contained in this Instructions for Use

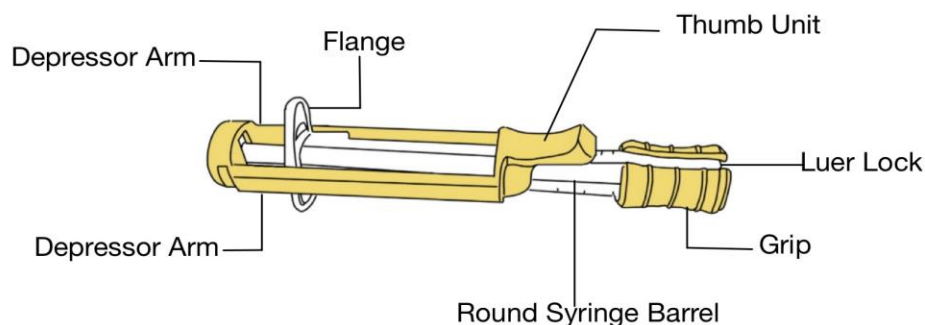

**Intended Use:** The Sales Syringe Adapter is intended to aid the user's dexterity in an applicable setting by enhancing the user's grip, control, and comfort in operating a syringe during one handed non-surgical procedures.

## Contents of Packaging

- Depressor
- Grip
- Instructions for use

## Contraindications

- Do not use this SSA for surgical procedures.

## Warnings ⚠

- Single use only.
- This device is to be used on a BD Luer-Lok 1ml/ 1 cc syringe only.
- Do not use in sterile procedures. As it may lead to Surgical Site Infection (SSI).

## How Supplied: 50

devices are packaged into a Snap Lock Box, clean, and non-sterile.

## Precautions


- Prior to use, inspect devices to ensure proper function and condition.
- Do not use the Syringe Adapter if it does not satisfactorily perform its intended function or if it has physical damage.
- Avoid mechanical shock or overstressing the devices.
- Caution needs to be stored at room temperature to avoid morphing of product.
- If cracked, sharp edges can harm user or patient.
- Caution, not designed to displace the full length of the syringe

# Sterile Sales Syringe Adapter

## Instructions for Use

### Directions for Use

1. To use a BD syringe Luer-Lok syringe, ensure all components are intact, functional, and packaged properly.
2. Position the barrel flange between the depressor arms and adjust the plunger for enhanced grip.
3. The "grip" component of the Sales Syringe Adapter should be snapped onto the barrel below the depressor arms.
4. Snap the barrel into the depressor by holding the round housing and syringe round plunger flange between the thumb and index fingers.

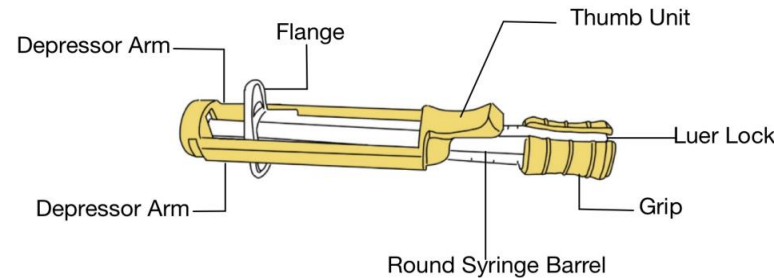

For a more detailed Instructions for Use:



MADE, The University of Iowa  
25 South Grand Ave.  
Iowa City, IA 52242  
(319) 384-3425

For product ordering or reporting of adverse events, please call (319) 384-3425, email [iowa-made@uiowa.edu](mailto:iowa-made@uiowa.edu), or visit [iowamade.org](http://iowamade.org)

Rev. B (09/01/2023)  
ID0003