

Sterile Sales Syringe Adapter

Instructions for Use

Important

- This product is intended to be used by a professional with appropriate training and experience
- Do not attempt to use the Sales Syringe Adapter (SSA) before completely reading and understanding the information contained in this Instructions for Use

Contents of Packaging

- Depressor
- Grip
- Instructions for use

How Supplied: Devices are individually packaged, clean, and sterile.

Contraindications


- Not intended for surgeries that consist of metered dosages.

Precautions

- ☐ Prior to use, inspect devices to ensure proper function and condition.
- ☐ Do not use the Syringe Adapter if it does not satisfactorily perform its intended function or if it has physical damage.
- ☐ Avoid mechanical shock or oversteering the devices.
- ☐ Caution, needs to be stored at room temperature to prevent deformation of device.
- ☐ Caution, not designed to displace the full length of the syringe
- ☐ If cracked, sharp edges can harm user or patient.

Warnings ⚠

- ☐ Single use only.
- ☐ Do not re-sterilize, as it could deform the device.
- ☐ This device is to be used on a BD Luer-Lok 1ml/1. cc syringe only. Incorrect type of syringe could result in dropping the syringe while injecting.
- ☐ If the "grip" is used, attach it around the Luer Lock region before attaching a hypodermic syringe needle. If done in contrary, can result in risk of needle puncture.
- ☐ The Sales Syringe Adapter is not a replacement for proper training and technique.
- ☐ In surgical procedures, do not use non-sterile device

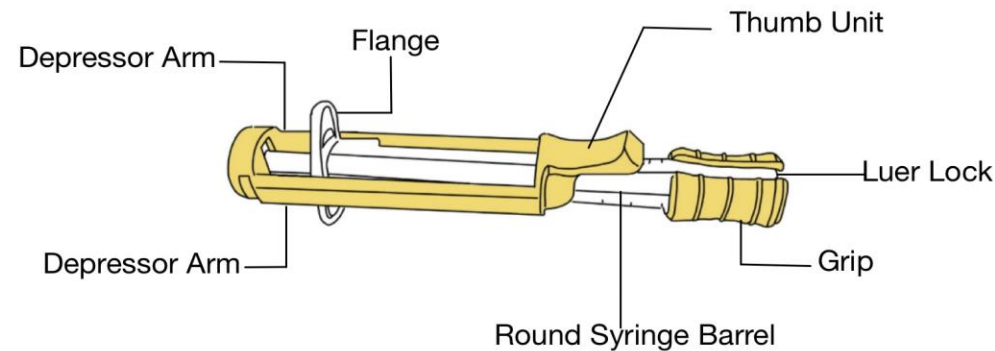


Intended Use: The Sales Syringe Adapter is intended to aid the user's dexterity in an applicable setting by enhancing the user's grip, control, and comfort in operating a syringe during one handed microsurgery procedures.

Sterile Sales Syringe Adapter

Instructions for Use

Directions for Use

1. To use a BD syringe Luer-Lok syringe, ensure all components are intact, functional, and packaged properly.
2. Position the barrel flange between the depressor arms and adjust the plunger for enhanced grip.
3. The "grip" component of the Sales Syringe Adapter should be snapped onto the barrel below the depressor arms.
4. Snap the barrel into the depressor by holding the round housing and syringe round plunger flange between the thumb and index fingers.

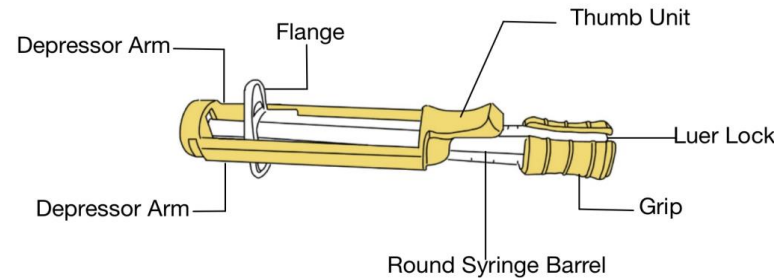

For a more detailed Instructions for Use:



MADE, The University of Iowa
25 South Grand Ave.
Iowa City, IA 52242
(319) 384-3425

For product ordering or reporting of adverse events, please call (319) 384-3425, email iowa-made@uiowa.edu, or visit iowamade.org

Rev. B (09/01/2023)
ID0003